

REPORT DOCUMENTATION PAGE				<i>Form Approved</i> OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Department of Defense, Washington Headquarters Services, Directorate for Information Operations and Reports (0704-0188), 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to any penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number. PLEASE DO NOT RETURN YOUR FORM TO THE ABOVE ADDRESS.					
1. REPORT DATE (DD-MM-YYYY) 20 May 2015		2. REPORT TYPE Journal Article		3. DATES COVERED (From – To) August 2005 – January 2015	
4. TITLE AND SUBTITLE Biopolymers for Sample Collection, Protection, and Preservation				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S) Iryna Sorokulova, Eric Olsen, Vitaly Vodyannoy				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) USAF School of Aerospace Medicine Public Health & Preventive Medicine Department/PHE 2510 Fifth St. Wright-Patterson AFB, OH 45433-7913				8. PERFORMING ORGANIZATION REPORT NUMBER AFRL-SA-WP-JA-2015-0036	
9. SPONSORING / MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSORING/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION / AVAILABILITY STATEMENT Distribution A: Approved for public release; distribution is unlimited. Case Number: 88ABW-2015-1850, 8 Apr 2015					
13. SUPPLEMENTARY NOTES					
14. ABSTRACT One of the principal challenges in the collection of biological samples from air, water, and soil matrices is that the target agents are not stable enough to be transferred from the collection point to the laboratory of choice without experiencing significant degradation and loss of viability. At present, there is no method to transport biological samples over considerable distances safely, efficiently, and cost-effectively without the use of ice or refrigeration. Current techniques of protection and preservation of biological materials have serious drawbacks. Many known techniques of preservation cause structural damages, so that biological materials lose their structural integrity and viability. We review applications of a novel bacterial preservation process, which is nontoxic and water soluble and allows for the storage of samples without refrigeration. The method is capable of protecting the biological sample from the effects of environment for extended periods of time and then allows for the easy release of these collected biological materials from the protective medium without structural or DNA damage. Strategies for sample collection, preservation, and shipment of bacterial, viral samples are described. The water-soluble polymer is used to immobilize the biological material by replacing the water molecules within the sample with molecules of the biopolymer. The cured polymer results in a solid protective film that is stable to many organic solvents, but quickly removed by the application of the water-based solution. The process of immobilization does not require the use of any additives, accelerators, or plastifiers and does not involve high temperature or radiation to promote polymerization.					
15. SUBJECT TERMS Biopolymer, acacia gum, trehalose, water, soil, air, pathogens					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT SAR	18. NUMBER OF PAGES 11	19a. NAME OF RESPONSIBLE PERSON Lt Col Eric Olsen
a. REPORT U	b. ABSTRACT U	c. THIS PAGE U			19b. TELEPHONE NUMBER (include area code)

Biopolymers for sample collection, protection, and preservation

Iryna Sorokulova¹ · Eric Olsen² · Vitaly Vodyanoy¹

Received: 8 April 2015 / Revised: 5 May 2015 / Accepted: 6 May 2015
© Springer-Verlag Berlin Heidelberg 2015

Abstract One of the principal challenges in the collection of biological samples from air, water, and soil matrices is that the target agents are not stable enough to be transferred from the collection point to the laboratory of choice without experiencing significant degradation and loss of viability. At present, there is no method to transport biological samples over considerable distances safely, efficiently, and cost-effectively without the use of ice or refrigeration. Current techniques of protection and preservation of biological materials have serious drawbacks. Many known techniques of preservation cause structural damages, so that biological materials lose their structural integrity and viability. We review applications of a novel bacterial preservation process, which is nontoxic and water soluble and allows for the storage of samples without refrigeration. The method is capable of protecting the biological sample from the effects of environment for extended periods of time and then allows for the easy release of these collected biological materials from the protective medium without structural or DNA damage. Strategies for sample collection, preservation, and shipment of bacterial, viral samples are described. The water-soluble polymer is used to immobilize the biological material by replacing the water molecules

within the sample with molecules of the biopolymer. The cured polymer results in a solid protective film that is stable to many organic solvents, but quickly removed by the application of the water-based solution. The process of immobilization does not require the use of any additives, accelerators, or plastifiers and does not involve high temperature or radiation to promote polymerization.

Keywords Biopolymer · Acacia gum · Trehalose · Water · Soil · Air · Pathogens

Introduction

Contaminated air, water, and soil bring the threat to health, life, and property. Efficient and timely sample collection that allows transferring agents from the environment to detectors or laboratories is vital for prevention and environmental remediation of the agent exposure. A fundamental knowledge of sample collection from various matrices is crucial.

Recovery and preservation of microorganisms prior to analysis are important considerations of overall detection and identification schemes. Collection of suspected pathogens from the environment is usually based on air filtration or swab sampling followed by microbial analysis using traditional culture or PCR assay (Temprano et al. 2004). However, at least one study shows that only 10–13 % of bacterial strains could be recovered using these methods (Meunier et al. 2005). The main problem is that microbiological analysis first relies on adequate on-hand storage methods or immediate transport of microbes in the laboratory in order to ensure viability (Rutala et al. 2006). Current techniques of recovery and preservation have serious drawbacks, including structural damage to cell walls and degradation of nucleic acids (Beal et al. 2001;

Electronic supplementary material The online version of this article (doi:10.1007/s00253-015-6681-3) contains supplementary material, which is available to authorized users.

✉ Vitaly Vodyanoy
vodyavi@auburn.edu

¹ Department of Anatomy, Physiology and Pharmacology,
Auburn, AL 36849, USA

² US Air Force School of Aerospace Medicine, Epidemiology
Laboratory, Wright-Patterson Air Force Base, Dayton, OH 45433,
USA

Broadbent and Lin 1999; Carvalho et al. 2004; Gouesbet et al. 2001). Consequences of freezing and freeze-drying include protein denaturation that leads to the decrease in cell viability (Carvalho et al. 2004; Leslie et al. 1995). The known processes of encapsulating immobilized bacteria in various starches involve spray-drying at elevated temperatures (Corcoran et al. 2004; Lian et al. 2002; O’Riordan et al. 2001) which is capable of damaging the microbes beyond repair. When preservation with conventional polymers is used, it usually involves elevating the temperature of the material or the use of UV radiation, both of which can again damage biological systems and nucleic acids. Resins are also used to some degree for the preservation of samples for electron microscopy, but the biological materials cannot be released for propagation after curing. Another method for encapsulating bacteria for use in biodegradation of gasoline involves a complex process using gellan gum (Moslemy et al. 2002). Many of these techniques have employed the use of protectants such as α - α trehalose, dimethyl sulfoxide (DMSO), glycerol, lactose, sucrose, and skim milk powder (Broadbent and Lin 1999; Carvalho et al. 2004; Conrad et al. 2000; Panoff et al. 2000), which contribute to the complexity of the preservation process. All of the aforementioned techniques are impractical for collecting samples in the field, are time-consuming, require special equipment and training, involve complex formulations, do not involve a broad spectrum of bacteria, and are limited in their applications.

Basic preservation methods, like refrigeration or tissue culture, have shortcomings such as limited shelf life, substantial cost, and possibility of contamination. As an alternative, cryopreservation has been effectively employed for long-term storage of biologicals. This method is founded on the concept that chemical, biological, and physical processes are adequately slowed down at cryogenic temperatures (-196°C) (Mazur 1984). Usually, cryopreservation protocols use large concentrations (e.g., 1.0–2.0 M) of cryoprotectants like DMSO, glycerol, or ethylene glycol (Mazur 1984) simply because they easily pass through cell membranes and so provide protection to the intracellular components. However, cooling media that contain DMSO have been demonstrated to be toxic (Abrahamsen et al. 2002). A substitute to DMSO consists of nontoxic sugars such as trehalose. Trehalose is a nontoxic, nonreducing disaccharide of glucose that has a capability to stabilize and preserve cells and cellular components throughout freezing (Hino et al. 1990; Sano et al. 1999) and drying (Chen et al. 2001). The application of trehalose as a cryoprotectant may be related to its protective interactions with lipid membranes and proteins during freezing-thawing (Carpenter et al. 1986). An impediment to the application of trehalose as a cryoprotectant is its impermeability to mammalian cell membranes (Chen

et al. 2001). Poration agents like ATP and benzoyl ATM need to be utilized to provide membrane permeability for trehalose (Buchanan et al. 2010).

Most of the existing methods of cell preservation are focused on the nucleic acid preservation and not on the safeguarding of cells (Gray et al. 2013). However, preservation of viable cells has definite advantages compared to the preservation of DNA, RNA only. First, the preserved cells could be multiplied by cultivation and provide large samples for testing, and second, live cells can also be used for physiological, pharmacological, and toxicological testing.

Preservation with biopolymers

Acacia gum has been used for thousands of years. In Ancient Egypt, it was used as a pigment binder and adhesive in paints for making hieroglyphs, binder in cosmetics and inks, and for adhering flaxen wrappings for embalming mummies (Sanchez et al. 2002). In current times, acacia gum is used extensively in food industry as an emulsifier, flavoring agent/adjuvant, formulation aid, stabilizer/thickener, humectant, surface-finishing agent; in pharmaceutical industry as a suspending agent, emulsifier, adhesive, and binder in tableting and in demulcent syrups; in cosmetics as stabilizer in lotions and protective creams, and an adhesive agent in blusher; foam stabilizer in liquid soaps; in lithography; in textile industry; in painting as dispersant in paints; and as dispersant in insecticidal/ acaricidal emulsions (Baldwin et al. 1999; Joseleau and Ullmann 1990; Sanchez et al. 2002). We discovered that one of the most important applications is the use of acacia gum in preserving microorganisms for extended periods of time without refrigeration (Krumnow et al. 2009; Sorokulova et al. 2008, 2012).

The preservation with acacia gum involves a single compound and simple procedure. The process of immobilization in acacia gum does not include any additives, accelerators, or plastifiers, nor does it involve elevated temperatures or radiation to promote polymerization. The polymerization and replacement of water occur spontaneously in a single process. The method is reversible, and the solid polymer can be safely removed from bacteria.

Acacia gum (or gum arabic) is an exudate collected from *Acacia senegal* stems and branches (Baldwin et al. 1999; Joseleau and Ullmann 1990). Most of the acacia gum (90 %) obtained for commercial use comes from Sudan, Africa (Joseleau and Ullmann 1990). The exudate is a hydrocolloid composed of various molecular components of D-galactose (~40 %), L-arabinose (~24 %), L-rhamnose (~13 %), D-glucuronic acid (~21 %), 4-O-methyl-D-glucuronic acid (~2 %), and polypeptide (~2 %) (Michel et al. 1998; Mocak et al. 1998).

Large aggregates have been fractionated which include an arabino-galactan protein ($2.3 \times 106 \text{ g mol}^{-1}$), arabino-galactan ($2.7 \times 105 \text{ g mol}^{-1}$), and glycoproteins (Baldwin et al. 1999; Sanchez et al. 2002). The structural unit for these aggregates is “ β -(1,3)-linked galactose backbone with branches of β -(1,6)-linked galactose containing arabinose, rhamnose, uronic acids and their derivatives” (Sanchez et al. 2002). Acacia gum readily dissolves in cold and hot water in concentrations up to 50 % (JECFA 1998). The solutions are characterized by a low viscosity and pH range 4.5–5.5. This polymer is indigestible to both humans and animals, not degraded in the intestine, but fermented in the colon by microflora (Michel et al. 1998). It is biodegradable and has no ecological impairments.

Another natural biopolymer, trehalose, is often used for the preservation of biologicals. Trehalose is a stable, non-reducing disaccharide with two linked glucose molecules. Trehalose exists naturally in plants, animals, and microorganisms and has long been consumed by humans as a component of mushrooms, baker's and brewer's yeasts, seaweeds and such invertebrates as lobsters. Trehalose plays a vital role in the preservation of biomembranes and the revival of certain biological functions following desiccation or freezing (Teramoto et al. 2008).

Comparison of different techniques used for the preservation of biologicals at ambient temperature is given in Table 1. The preservation with a SampleMatrix[®] offered by Biomatrix[®], San Diego, CA, is based on the proprietary synthetic compounds that mimic the protective properties of disaccharides, like trehalose (Lee et al. 2014). It is suitable for the preservation of DNA, RNA, and *Escherichia coli* strains only. It is not recommended for use with other bacterial species. DNAgard, another product of Biomatrix, was created to permeate the cells and preserve the DNA in a proprietary solution. This system does not require refrigeration or freezing to retain the stability of the DNA for long-term storage (Table 1) (Gray et al. 2013).

Recently, preservation of stem cells at room temperature for 4 weeks was demonstrated with trehalose (Buchanan et al. 2010). The method is very labor intensive and expensive. Beside the step of cell poration, this technique also includes a step of freeze-drying, which is not readily available in field conditions.

FTA technology (Whatman, GE Healthcare, Little Chalfont, UK) is based on the proprietary FTA cards that contain chemicals that lyse cells, denature proteins, and protect nucleic acids from damage (Table 1) (Gray et al. 2013). A solution-based mix of high salts, RNAlater[®] (Ambion, Carlsbad, CA) (Table 1), preserves tissues for a long time depending on the temperature. For extended storage, fresh samples are immersed in the solution, refrigerated overnight allowing saturation of the tissue, and frozen after that. RNAlater[®] has been demonstrated to

Table 1 Comparison of cell and nucleic acid preservation techniques

System	Medium	Preservation	Room temperature stability	Long-term storage	Reference
Acacia gum biopolymer	Biopolymer, dry	Gram-negative, Gram-positive bacteria, DNA, phages, antibodies, and proteins	Months, years	Room temperature	(Krumnow et al. 2009; Sorokulova et al. 2008, 2012)
Biomatrix [®]	SampleMatrix [®]	<i>E. coli</i> strains only. Not recommended for use with other bacterial species.	1 month for viable <i>E. coli</i>	Room temperature	(Biomatrix 2013)
Stabilizing sugar	Trehalose, liquid	DNA, RNA	4 weeks	Room temperature	(Buchanan et al. 2010)
FTA [®]	Filter card	Stem cells	Years	Room temperature	(Gray et al. 2013)
FTA elute [®]	Filter card	DNA	Years	Room temperature	(Gray et al. 2013)
DESS	Liquid	DNA	Months	Room temperature	(Gray et al. 2013)
RNAlater [®]	Liquid	RNA	1 week 1 month Years	Room temperature, 4 °C–20 °C	(Gray et al. 2013)
DNAgard™	Liquid	DNA	1 month	Room temperature	(Gray et al. 2013)

provide greater DNA yield than FTA® cards. RNAlater® furthermore gives the efficient preservation of RNA (Gray et al. 2013). DMSO/EDTA/saturated sodium chloride (DESS) is a nonproprietary solution that has been used to preserve DNA.

It was found that the liquid-based preservatives (DNAgard™, RNAlater®, and DESS) performed better than the card-based methods. Not one liquid technique evidently outperformed the others (Gray et al. 2013). An acacia gum biopolymer shows a very broad spectrum of preservation, including Gram-negative, Gram-positive bacteria, DNA, phages, antibodies, and proteins. It is also simple compared to other methods.

Among the methods available for assessing airborne viable particles (sieve/nozzle impactor, slit-to-agar, centrifugal, filtration, and impinge), the filtration may be the most suitable method for microbiological monitoring (Temprano et al. 2004). Specialized filters are the most suitable for the collection and analysis of airborne agents. Similarly, agents from water and soil are to be processed and filtered (Roger 1994).

Analysis of literature data allows us to conclude that method of preservation with acacia gum-based biopolymer is much more superior compared to other known to proposer's methods. It is simple, efficient, and cost-effective.

Examples of preservation of microorganisms in natural biopolymers

Researchers used different methods to characterize the rate of the cell population reduction with time.

In 1943, (Katzin et al. 1943) defined the decimal reduction (D-value) time as follows:

$$D = \frac{2.3}{k} = t / \log \left(\frac{C_1}{C_2} \right) \quad (1)$$

where k is the monomolecular reaction rate constant, and C_1 and C_2 are the initial and final concentrations of bacteria, respectively, subjected to a constant lethal temperature for $t=(t_2-t_1)$ duration of time. In this model, the classical D-value presents a simple biological significance: time that leads to a tenfold reduction of surviving population and is easily estimated from a simple linear regression. Other population reduction values can be readily converted to the D-value for comparison. The bacterial survival ratio (BSR) defined as $100 \times \log C_1 / \log C_2$ (Munoz-Rojas et al. 2006) can be transformed into the D-value of the following equation:

$$D = t_{\text{BSR}} / \log C_1 (1 - \text{BSR} / 100) \quad (2)$$

The D-value is related to the time to degrade to 100 cells (Sorokulova et al. 2012), t_{100} by the following expression:

$$D = t_{100} / (\log C_1 - 2) \quad (3)$$

The D-value calculated from the percent of survival $S=(C_2/C_1)\%$

$$D = t / \log \left(\frac{100}{S\%} \right) \quad (4)$$

The natural biopolymer, trehalose, is frequently utilized for the preservation of microorganisms. Trehalose significantly enhanced the tolerance of *E. coli* to air-drying (Table 1). The protective effects of trehalose during desiccation appear to be due to its stabilizing effect on membrane structure, its chemically inert nature, and the tendency of trehalose solutions to form glasses upon drying (Welsh and Herbert 1999).

The stability of a vacuum-dried *Lactobacillus paracasei* F19 was characterized by drying process in the presence of trehalose (Foerst et al. 2012). The D-value of the stored bacteria with trehalose was about four times larger compared with a D-value of unprotected cells (Table 2).

The protection capabilities of dry biopolymers were illustrated by the accelerated aging in alginate, POLYOX, carrageenan, and acacia gum (Rojas-Tapias et al. 2013, 2015). They reported storage enhancement of 1.8–5.5 (Table 2).

Salmonella typhimurium was immobilized in the solid polymer. After 21 days at room temperature, the bacteria were released and showed no difference in motility and viability compared to initial culture (Figure S1, and Preservation video, Supplementary Materials).

The formulations with *Bacillus subtilis* spores were tested for spore stability during storage at temperatures ranging from 40 to 90 °C and for bacterial release. Thermodynamic analysis showed that immobilization of spores in acacia gum significantly increased their viability compared with that of unprotected spores. The viability was further enhanced when suspensions of spores in acacia gum were added to the granules of charcoal (Table 2, Fig. 1). The number of the spores released after storage was also increased when spores were treated with acacia gum prior to immobilization in tapioca and charcoal. Formulations of *Bacillus* spores with acacia gum and porous carriers (charcoal and tapioca) prolong the anticipated shelf life of spores even under ambient temperature and provide slow and steady bacterial release consistent with their high viability (Sorokulova et al. 2008).

We demonstrated and characterize the efficacy of acacia gum biopolymer for preservation of pathogenic bacteria (*Bacillus anthracis* and methicillin-resistant *Staphylococcus aureus* (MRSA)) on different materials, used for swabbing and filtration: cotton, wool, polyester, rayon, charcoal cloth, and Whatman paper. The biopolymer used for the preservation of two pathogens has been shown to significantly

Table 2 Preservation and storage of microorganisms immobilized in dry biopolymers

Biopolymer	Microorganism	Temperature °C	D-value ^a	Reference
Trehalose	<i>E. coli</i> NCIB 9484, unprotected/protected	25	3.3/36.3 days	(Welsh and Herbert 1999)
Trehalose	<i>Lactobacillus Paracasei</i> Unprotected/protected	37	6.6/25.5 days	(Foerst et al. 2012)
Alginate	<i>Azotobacter chroococcum</i> Unprotected/protected	15	6/23.7 days	(Rojas-Tapias et al. 2013, 2015)
POLYOX ^{®b}	<i>Azotobacter chroococcum</i>	15	24.5 days	(Rojas-Tapias et al. 2013)
Carrageenan	<i>Azotobacter chroococcum</i>	15	32.9 days	(Rojas-Tapias et al. 2013)
Acacia gum	<i>Azotobacter chroococcum</i>	15	10.6 days	(Rojas-Tapias et al. 2013)
Acacia gum	<i>E. coli</i> O1:K1:H7 (ATCC 11775) Unprotected/protected	15	8/57.5 days	(Krumnow et al. 2009)
Pullulan	<i>E. coli</i> O1:K1:H7 (ATCC 11775) Unprotected/protected	15	8/33 days	(Krumnow et al. 2009)
Pullulan	<i>Bacillus subtilis</i> , ATCC 6051, humidity 76 %	25	52.3 days	(Krumnow et al. 2009)
Acacia gum	<i>Bacillus subtilis</i> , ATCC 6051, humidity 76 %	25	272 days	(Krumnow et al. 2009)
Acacia gum	<i>Bacillus subtilis</i> , ATCC 6051. Unprotected spores/Spores immobilized in charcoal+acacia gum	20	1.6/48 years	(Sorokulova et al. 2008)
Acacia gum	<i>S. aureus</i> resistant to methicillin (MRSA) immobilized in filter paper/Filter paper+acacia gum	15	43/89 days	(Sorokulova et al. 2012)
Acacia gum	<i>B. anthracis</i> Sterne spores, immobilized in filter paper/Filter paper+acacia gum	50	1.1/3.1 years	(Sorokulova et al. 2012)

^a D-value is a time that leads to a tenfold reduction of surviving population

^b POLYOX, water-soluble resins from Colorcon (Harleysville, USA)

protect bacteria during dehydration and storage in all tested samples at the range of temperatures (5–45 °C for MRSA and 40–90 °C for *B. anthracis*). Our results showed higher recovery as well as higher viability during the storage of both bacteria in all materials with acacia gum (Table 2). The addition of acacia gum polymer to swabbing materials or filters will increase the efficacy of sample collection and identification of pathogenic bacteria from locations such as hospitals or the environment. The approach can also be used for long-term storage of culture collections since acacia gum contributes to the viability and stability of bacterial cultures (Sorokulova et al. 2012).

Natural polymers of acacia gum and pullulan were used to preserve model bacteria *E. coli* and *B. subtilis* via immobilization and storage under various conditions. Formulation of *E. coli* and *B. subtilis* in acacia gum significantly increased the viability of both cultures during desiccation at 40 °C as well as during the storage at various temperatures and relative humidity. In the ranges of temperatures and humidity used in experiments, the high humidity affected the viability of *E. coli* more than high temperature. Thermodynamic parameters for *E. coli* thermal degradation were used for quantification of the results and characterization of the preservation process. The viability of *B. subtilis* in the acacia gum polymer was not significantly

changed during the storage in the temperature and humidity experiments. The number of viable *B. subtilis* recovered after storage in pullulan, and in PBS under various humidity conditions, was 1–2 logs less in comparison with the number of cells before storage. It was found that acacia gum provides better protection than pullulan for both bacteria during the preservation process (Table 2) (Krumnow et al. 2009).

Sampling and preservation strategies

Standard and enhanced procedures for sampling collection

Stability of microbial cultures during sampling and storage is a vital issue in various fields of medicine, biotechnology, food science, and forensics. Acacia gum biopolymer-based preservation has many advantages. It is a simple and inexpensive process that does not require specialized equipment or training to preserve bacteria successfully. The most important application of the biopolymer is collecting samples in the field. It has been shown in our preliminary studies that acacia gum is capable of preserving bacteria in stringent conditions such as high temperature and humidity without refrigeration. Thus,

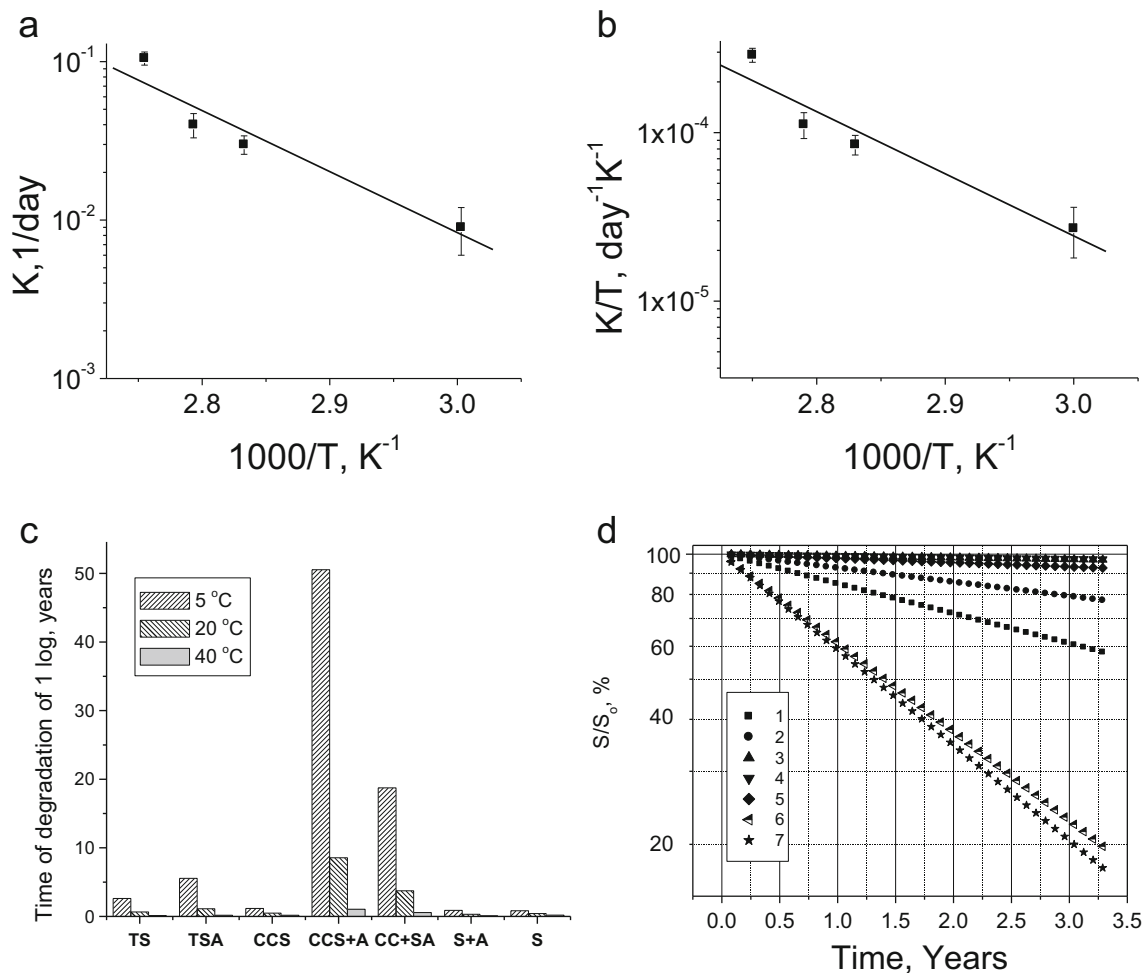


Fig. 1 The estimated viability of spores by Arrhenius and Eyring equations. **a** Points are experimental data plotted by Arrhenius plot; the line is a linear regressions: $R(P) = -0.96$ (<0.04). **b** Points are experimental data plotted by Eyring equation, and the line is a linear regressions. $R(P) = -0.96$ (<0.05). **c** Time for degradation of one log of spores at different temperatures. **d** The viability of spores (percent of

surviving spores, S/S_0 %) at 5 °C; 1—tapioca+spores (TS); 2—tapioca+(spores+acacia gum) (T+SA); 3—charcoal+spores (CCS); 4—(charcoal+spores)+acacia gum (CCS+A); 5—charcoal+(spores+acacia gum) (CC+SA); 6—spores+acacia gum (SA); 7—spores (S) (from Sorokulova et al. 2008, by permission from WILEY # 3603090702078)

collecting samples at remote locations using biopolymer may be applied to preserving microorganisms for further testing in a laboratory. Normally, the smaller the time that elapses between the collection of a sample and its examination, the more reliable will be the analytical results. The Table 1 shows that acacia gum biopolymer is the most effective compared to other cell preservation methods. Additionally, the acacia gum biopolymer is the only one that was examined to preserve a large variety of biological.

The standard procedure followed the environmental sample collection that includes first refrigeration by ice packs, dry ice, or portable refrigerator, then transportation and storage, and finally the laboratory analysis (Kelly et al. 2006; NEWWA 2008; ODH 2008). The enhanced method of the sample collection replaces the step of refrigeration with the preservation utilizing the acacia gum biopolymer. The enhancement of the modified method can be illustrated by a

simple example. If one collects 1000 water samples according to a general First Responder Procedure, water tubes with samples are cooled to +4 °C with ice or refrigeration and transferred to the laboratory for storage and testing. The volume of the sampling containers may be 40 mL, 60 mL, or 1 L, depending on the analytical method used (NEWWA 2008). Cooler boxes need to be filled with ice to have ice to sample ratio of at least 2:1 to achieve an ice retention for about 3–4 days (YETI 2013). The price of a single 60-mL sample is a sum of the costs of a bottle, ice, freight transportation, and storage. Five percent of a container weight needs to be added to a weight of a sample. Using acacia gum biopolymer method, a nonwater content of 60 mL sample is filtered. The filter with microorganisms is immobilized in biopolymer and transferred to a test location in a small zip bag. Table 3 compares the general method and the method that utilizes the acacia gum biopolymer for preservation during the transport and

Table 3 Comparison of sample collection price and efficiency by different methods

Properties	Methods		Gain by the proposed method
	Standard	Acacia gum biopolymer	
Disposables ^a , \$ per sample	0.72	0.78	0.92
Ice ^b /biopolymer ^c , \$ per sample	0.048	0.01	4.8
US Freight by truck ^d , 1000 miles, \$ per sample	0.015	0.000074	200
Efficiency ^e , \$ per sample per day	0.0979	0.0376	2.6

^a 60 mL bottles (Dynalox 2011) are used in the standard method. In the biopolymer method, 2×3 zipper seal bags (Dynalox 2011), filters (VWR #28157-949), and Petri dishes (VWR #82050-536)

^b \$0.4 is a price of 1 kg of ice. One sample requires 120 g of ice, and it costs \$0.048

^c The cost of biopolymer in one sample is \$0.01 (Sigma-Aldrich #G9752)

^d It takes 3357 BTU to move 1 t 1000 miles (Davis et al. 2014), and it costs \$21.98 per million BTU (EIA 2012). The weight of one cooled sample in the standard method (W_1) is comprised of the weight of the empty bottle, sample, ice, and container. $W_1 = 13 + 60 + 120 + 9.7 \approx 204$ g. The sample weight in the biopolymer method (W_2) is comprised of the weight of the empty bag, filter, solid biopolymer, ice, and container. $W_2 = 0.4 + 0.1 + 0.2 + 0.35 \approx 1$ g. Therefore, US Freight by truck, 1 mile for the standard method costs = $3357 \times 204 \times \$21.98 \times 1.10^{-9} = \0.015 , while in the proposed method costs = $3357 \times 1 \times \$21.98 \times 1.10^{-9} \approx \0.000074 .

^e The efficiency of sample collection takes into account a viability loss during sampling, transfer, and storage. The higher the number of viable microorganisms, the greater efficiency of sampling can be. The efficiency is exemplified by a sampling of *E. coli*. At 15 °C, D-values with and without biopolymer are 57 and 8 days, respectively, (Table 2). The efficiency is calculated as total price per sample/D-value

storage before tests. The sampling efficiency of particular microorganism is defined by the price of sampling per sample per D-value of the microorganism.

As can be seen from the Table 3, this preservation system significantly improves the performance of the environmental sampling. The collected samples need no refrigeration for transport and storage. The weight of collected samples is dramatically reduced. One thousand 60 mL samples on ice weigh 204 kg, while the same number samples preserved by biopolymer weigh only 1 kg. Therefore, the freight is about 200 less expensive. The efficiency of sampling (the price of collecting and transfer viable bacteria) is strongly increased. As exemplified by sampling of *E. coli* in water, the sampling efficiency with biopolymer is increased by 260 %. The use of filters as a carrier for biopolymer prolongs the viability of microorganisms at the average by four times (Krumnow et al. 2009). Therefore, the efficiency of sampling with carrier filters and biopolymer may over-perform the standard method by more than ten times.

Concept of preservation at ambient temperature by biopolymer

Acacia gum biopolymer is a natural product that can preserve cells and other biologicals for extended periods at room temperature. Current methods require cold temperatures, which are labor intensive, expensive to maintain, and make the samples difficult to transport. The biopolymer preserves embed samples in a stable, protective film. Samples can be recovered by simply adding water, without the need for heat or radiation. Samples have been kept alive for over 6 months at room

temperature and over 16 months at refrigerated temperatures (Krumnow et al. 2009; Sorokulova et al. 2008, 2012)

The biopolymer technology is applied in three stages:

1. Immobilization and polymerization
 - Biopolymer provides structural integrity, guards against mechanical stress, and slows metabolic processes in cells.
2. Stabilized state existence (dormant state for cells)
 - The dormant state is characterized by slow metabolism, water replacement; partial water retention; and molecular stability.
3. Release
 - Biopolymer protects against structural changes, releases mechanical stress, and facilitates the increase of metabolism in cells.

Microorganism's collection and preservation schemes for contaminated air, water, and soil

Air samples

The samples of the environmental air are collected by filtering relatively large volumes of air drove by a battery-driven

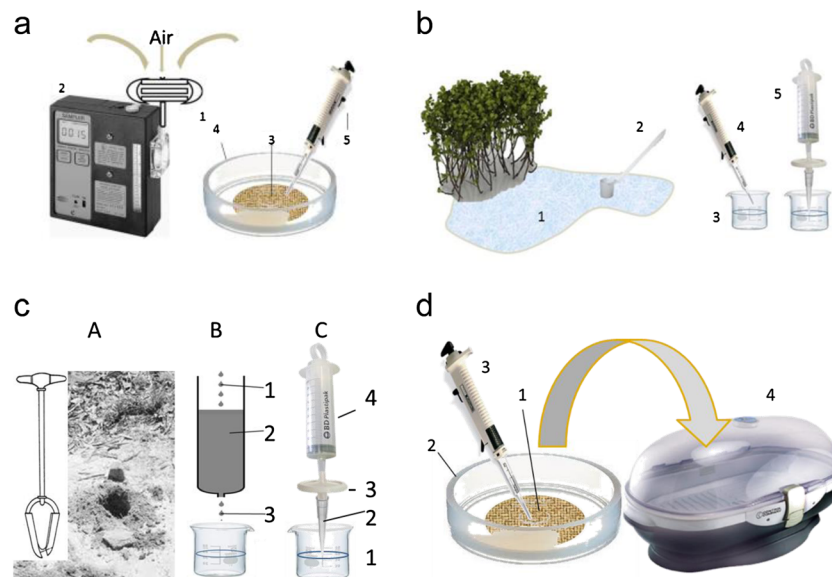


Fig. 2 Sample collection from the environment. **a** Sample collection of environmental air. 1—filter cassette with filter, 2—portable battery-driven pump, 3—filter with collected bacteria from ambient air, 4—Petri dish, 5—micropipette with liquid biopolymer. **b** Sample collection of environmental water. 1—body of water, 2—dipper, 3—beaker with collected environmental water, 4—syringe filtration system. **c** Collection of environmental soil samples. **A** represents the collection of soil samples

with an auger. **B** represents the microorganism elution (1—water, 2—soil, 3—suspension of microorganisms). **C** represents the collection of microorganisms by syringe filtration. **d** Preservation of microorganisms with biopolymer. 1—filter with collected microorganisms, 2—Petri dish, 3—micropipette with liquid biopolymer, 4—battery-driven potable desiccator

portable pump as shown in Fig. 2a. The filter with trapped environmental microorganisms (3) is positioned at the bottom of the sterilized Petri dish (4) and added liquid biopolymer. After the last step, the sample is ready for polymerization.

Water samples

The environmental samples of water are collected from a natural body of water (1) using a dipper (2) as shown in Fig. 2b. The sample is added to the volumetric beaker (3) and filtered through a syringe filter system (4). After filtration, the filter with collected microorganisms is covered with liquid biopolymer.

Soil samples

The soil samples are collected from the environmental soil (Fig. 2c A). The samples then are inoculated by the sterile dH₂O (Fig. 2c B) and, finally, subjected to filtration and collection by the filter (Fig. 2c C). After that, the filter with the agent will be immobilized in biopolymer.

Sample polymerization

After microorganisms are collected on a filter, and each filter is positioned in the individual small Petri dishes, covered with liquid biopolymer, the Petri dishes with samples then are

covered with vented covers and transferred to the portable battery-driven desiccator to dry and polymerize for ~2 h. After that, the samples are ready for storage and analytical tests. The control filters with microorganisms are also going through drying in a desiccator (Fig. 2d).

Release of microorganisms and DNA analysis

After samples arrive at the laboratory, microorganisms are released and analyzed. Every sample can be placed in 10 mL of water buffer solution, homogenized, and filtered. The filtrate containing microorganisms is subjected to DNA analysis (Ogram et al. 1987; Park and Crowley 2005; Zhou et al. 1996).

Conclusions

The method of collection and preservation of microorganisms with the use of water-soluble biopolymers has a few advantages as compared with standard methods utilizing refrigeration.

- Immobilization of microorganisms in dry biopolymer allows the collection of biological samples from air, water, and soil matrices and transfers them from the collection

point to the laboratory of choice without experiencing significant degradation and loss of viability.

- Dry biopolymers are capable of preserving biological samples for extended periods (30+ days) without refrigeration.
- Water-soluble biopolymers ensure the complete recovery of preserved biological materials.
- The method of preservation with dry biopolymers is simple and rapid.
- Collected biological specimens are released from the biopolymer by the water-based solution.
- Hardens/solidifies under simple conditions (no temperature above 40 °C, no chemicals or radiation).
- Nontoxic and environmentally friendly.
- Economically sound.
- Acacia gum does not interfere with established sample examination techniques (microscopy, PCR, pulsed-field electrophoresis, and other methods.)
- Formulations of microorganisms with biopolymer and porous carriers (charcoal and tapioca) prolong the shelf life under ambient temperature and provide slow and steady release.
- Filters made of cotton, wool, polyester, rayon, charcoal cloth, and Whatman paper impregnated with acacia gum were shown to enhance the preservation of pathogenic bacteria actively.

Acknowledgments This study was supported by Air Force CRADA grant 13-268-SG-C13016. The views expressed in this article are those of the authors and do not reflect the official policy or position of the United States Air Force, Department of Defense, or the US Government.

Conflict of interest The authors declare that they have no competing interests.

References

- Abrahamsen JF, Bakken AM, Bruserud O (2002) Cryopreserving human peripheral blood progenitor cells with 5-percent rather than 10-percent DMSO results in less apoptosis and necrosis in CD34+ cells. *Transfusion* 42(12):1573–1580
- Baldwin TC, Quah PE, Menzies AR (1999) A serotaxonomic study of Acacia gum exudates. *Phytochemistry* 50(4):599–606
- Beal C, Fonseca F, Corrieu G (2001) Resistance to freezing and frozen storage of *Streptococcus thermophilus* is related to membrane fatty acid composition. *J Dairy Sci* 84(11):2347–2356
- Biomatrix (2013) Biostability <http://www.biomatrix.com/index.html>. Publisher
- Broadbent JR, Lin C (1999) Effect of heat shock or cold shock treatment on the resistance of *Lactococcus lactis* to freezing and lyophilization. *Cryobiology* 39(1):88–102
- Buchanan SS, Pyatt DW, Carpenter JF (2010) Preservation of differentiation and clonogenic potential of human hematopoietic stem and progenitor cells during lyophilization and ambient storage. *PLoS ONE* 5(9):e12518. doi:10.1371/journal.pone.0012518
- Carpenter JF, Hand SC, Crowe LM, Crowe JH (1986) Cryoprotection of phosphofructokinase with organic solutes: characterization of enhanced protection in the presence of divalent cations. *Arch Biochem Biophys* 250(2):505–512
- Carvalho AS, Silva J, Ho P, Teixeira P, Malcata FFXF, Gibbs P (2004) Relevant factors for the preparation of freeze-dried lactic acid bacteria. *International Dairy Journal* In Press, Corrected Proof
- Chen T, Acker JP, Eroglu A, Cheley S, Bayley H, Fowler A, Toner M (2001) Beneficial effect of intracellular trehalose on the membrane integrity of dried mammalian cells. *Cryobiology* 43(2):168–181. doi:10.1006/cryo.2001.2360
- Conrad PB, Miller DP, Cielenski PR, de Pablo JJ (2000) Stabilization and preservation of *Lactobacillus acidophilus* in saccharide matrices. *Cryobiology* 41(1):17–24
- Corcoran BM, Ross RP, Fitzgerald GF, Stanton C (2004) Comparative survival of probiotic *Lactobacilli* spray-dried in the presence of prebiotic substances. *J Appl Microbiol* 96(5):1024–1039
- Davis SC, Diegel SW, Boundy RG (2014) Transportation energy data book: Edition 33 U.S. Department of Energy cta.oeml.gov/data
- Dynalab (2011) Labware http://www.dynalabcorp.com/documents/Dynalab2011_LR.pdf
- EIA (2012) Annual energy review <http://www.eia.gov/totalenergy/data/annual/showtext.cfm?t=ptb0303>. US Energy Information Administration
- Foerst P, Kulozik U, Schmitt M, Bauer S, Santivarangkna C (2012) Storage stability of vacuum-dried probiotic bacterium *Lactobacillus paracasei* F19. *Food Bioprod Process* 90(2):295–300. doi:10.1016/j.fbp.2011.06.004
- Gouesbet G, Jan G, Boyaval P (2001) *Lactobacillus delbrueckii ssp. bulgaricus* thermotolerance. *Lait* 81:301–309
- Gray MA, Pratte ZA, Kellogg CA (2013) Comparison of DNA preservation methods for environmental bacterial community samples. *FEMS Microbiol Ecol* 83(2):468–477. doi:10.1111/1574-6941.12008
- Hino A, Mihara K, Nakashima K, Takano H (1990) Trehalose levels and survival ratio of freeze-tolerant versus freeze-sensitive yeasts. *Appl Environ Microbiol* 56(5):1386–1391
- JECFA (1998) Gum arabic. Food and nutritional paper 52(Addendum 6). FAO, Rome
- Joseleau J-P, Ullmann G (1990) Biochemical evidence for the site of formation of gum arabic in *Acacia senegal*. *Phytochemistry* 29(11):3401–3405
- Katzin LI, Sandholzer LA, Strong ME (1943) Application of the decimal reduction time principle to a study of the resistance of *Coliform* bacteria to pasteurization. *J Bacteriol* 45(3):265–272
- Kelly CD, Egan C, Cirino NM (2006) The CODE RED solution: bioterror response training for first responders. *Biosecure Bioterror* 4(4):391–396. doi:10.1089/bsp.2006.4.391
- Krumnow AA, Sorokulova IB, Olsen E, Globa L, Barbaree JM, Vodyanoy VJ (2009) Preservation of bacteria in natural polymers. *J Microbiol Methods* 78(2):189–194
- Lee SB, Crouse CA, Lline MC (2014) Optimizing storage and handling of DNA extracts. In: Shewale JG, Liu EH (eds) *Forensic DNA analysis*. CRC Press, Boca Raton, pp 19–38
- Leslie S, Israeli E, Lighthart B, Crowe J, Crowe L (1995) Trehalose and sucrose protect both membranes and proteins in intact bacteria during drying. *Appl Environ Microbiol* 61(10):3592–3597
- Lian WC, Hsiao HC, Chou CC (2002) Survival of bifidobacteria after spray-drying. *Int J Food Microbiol* 74(1–2):79–86
- Mazur P (1984) Freezing of living cells: mechanisms and implications. *Am J Physiol* 247(3 Pt 1):C125–C142
- Meunier O, Hernandez C, Piroird M, Heilig R, Steinbach D, Frey A (2005) Bacteriological sampling of surfaces: importance of the step of enrichment and the choice of the culture media. *Ann Biol Clin (Paris)* 63(5):481–486

- Michel C, Kravtchenko TP, David A, Gueneau S, Kozlowski F, Cherbut C (1998) In vitro prebiotic effects of Acacia gums onto the human intestinal microbiota depends on both botanical origin and environmental pH. *Anaerobe* 4(6):257–266
- Mocak J, Jurasek P, Phillips GO, Varga S, Casadei E, Chikemai BN (1998) The classification of natural gums. X. Chemometric characterization of exudate gums that conform to the revised specification of the gum arabic for food use, and the identification of adulterants. *Food Hydrocoll* 12(2):141–150
- Moslemy P, Neufeld RJ, Guiot SR (2002) Biodegradation of gasoline by gellan gum-encapsulated bacterial cells. *Biotechnol Bioeng* 80(2):175–184
- Munoz-Rojas J, Bernal P, Duque E, Godoy P, Segura A, Ramos JL (2006) Involvement of cyclopropane fatty acids in the response of *Pseudomonas putida* KT2440 to freeze-drying. *Appl Environ Microbiol* 72(1):472–477. doi:10.1128/aem.72.1.472-477.2006
- NEWWA (2008) Sampling guide for first responders to drinking water contamination threats and incidents <http://www.epa.gov/region1/eco/drinkwater/pdfs/1stRespondersGuide.pdf>. New England Water Works Association
- O’Riordan K, Andrews D, Buckle K, Conway P (2001) Evaluation of microencapsulation of a *Bifidobacterium* strain with starch as an approach to prolonging viability during storage. *J Appl Microbiol* 91(6):1059–1066
- ODH (2008) Environmental sample collection guide http://www.wvdhhr.org/labservices/shared/docs/bt/oh_state_laboratory_submission_guidelines.pdf. Ohio department of health laboratories
- Ogram A, Sayler GS, Barkay T (1987) The extraction and purification of microbial DNA from sediments. *J Microbiol Methods* 7(2–3):57–66. doi:10.1016/0167-7012(87)90025-X
- Panoff J-M, Thammavongs B, Gueguen M (2000) Cryoprotectants lead to phenotypic adaptation to freeze-thaw stress in *Lactobacillus delbrueckii ssp. bulgaricus* CIP 101027T. *Cryobiology* 40(3):264–269
- Park JW, Crowley DE (2005) Normalization of soil DNA extraction for accurate quantification of target genes by real-time PCR and DGGE. *Biotechniques* 38(4):579–586
- Roger S (1994) Guide to sampling air, water, soil, and vegetation for chemical analysis. In: State of California Environmental Protection Agency DoPR (ed). vol EH 94–04. Environmental hazards Assessment Program, Sacramento, California, p 1–57
- Rojas-Tapias D, Ortiz-Vera M, Rivera D, Kloepper J, Bonilla R (2013) Evaluation of three methods for preservation of *Azotobacter chroococcum* and *Azotobacter vinelandii*. *Univ Sci* 18:129–139
- Rojas-Tapias D, Sierra OO, Botía DR, Bonilla R (2015) Preservation of *Azotobacter chroococcum* vegetative cells in dry polymers. *Univ Sci* 20(2):201–207
- Rutala WA, White MS, Gergen MF, Weber DJ (2006) Bacterial contamination of keyboards: efficacy and functional impact of disinfectants. *Infect Control Hosp Epidemiol* 27(4):372–377
- Sanchez C, Renard D, Robert P, Schmitt C, Lefebvre J (2002) Structure and rheological properties of acacia gum dispersions. *Food Hydrocoll* 16(3):257–267
- Sano F, Asakawa N, Inoue Y, Sakurai M (1999) A dual role for intracellular trehalose in the resistance of yeast cells to water stress. *Cryobiology* 39(1):80–87. doi:10.1006/cryo.1999.2188
- Sorokulova I, Watt J, Olsen E, Globa L, Moore T, Barbaree J, Vodyanoy V (2012) Natural biopolymer for preservation of microorganisms during sampling and storage. *J Microbiol Methods* 88:140–146
- Sorokulova IB, Krumnow AA, Pathirana S, Mandell AJ, Vodyanoy V (2008) Novel methods for storage stability and release of *Bacillus* spores. *Biotechnol Prog* 24(5):1147–1153
- Temprano G, Garrido D, Daquino M (2004) Comparative study of airborne viable particles assessment methods in microbiological environmental monitoring. *PDA J Pharm Sci Technol* 58(4):215–221
- Teramoto N, Sachinvala ND, Shibata M (2008) Trehalose and trehalose-based polymers for environmentally benign, biocompatible and bioactive materials. *Molecules* 13(8):1773–1816
- Welsh DT, Herbert RA (1999) Osmotically induced intracellular trehalose, but not glycine betaine accumulation promotes desiccation tolerance in *Escherichia coli*. *FEMS Microbiol Lett* 174(1):57–63
- YETI (2013) Ice retention <http://www.yeticoolers.com/pages/FAQ.html>.
- Zhou J, Bruns MA, Tiedje JM (1996) DNA recovery from soils of diverse composition. *Appl Environ Microbiol* 62(2):316–322